

# URIGLOW® ORDERING INFORMATION

## Transilluminating Ureteric Stents & Light Guide Coupler

URIGLOW® Transilluminating Ureteric Stents provide rapid identification of the lower pelvic ureter during open and laparoscopic surgery.

**HIGH VISIBILITY** in a wide range of operative conditions has been achieved by using a new type of optical fibre, ensuring better intra-operative performance and easier detection.

**THE URIGLOW® LIGHT GUIDE COUPLER (LGC)** is a reusable precision optical device designed to absorb up to 90% of infrared radiation present in the output of medical light sources. The URIGLOW® LGC safely and securely links fibre optic cables to the URIGLOW® Transilluminating Ureteric Stents.

**HIGH LEVEL IR ABSORPTION** The URIGLOW® LGC is designed to absorb up to 90% of light source IR output which means there is no detrimental heating effect in the ureter from high intensity light sources.

**ROBUST DESIGN** Solid stainless steel construction securely protects the delicate IR mirror and lens assembly.

**SECURE ATTACHMENT** The URIGLOW® LGC will securely attach the URIGLOW® Transilluminating Ureteric Stents to any standard Storz/ACMI screw type light cable fitting.

**ROBUST DESIGN** The URIGLOW® LGC is designed for steam sterilisation. Autoclave; 134°C (+3°C - 0°C) for a minimum of three minutes. Carefully follow instructions provided with the unit.



### SPECIFICATIONS

#### Optical fibre:


1.9mm (6Fr) OD x 100cm radiopaque marker line.

#### Active tip:

6 x 1cm high intensity emission points.  
1st point 15mm from distal domed tip.

#### Marker positions from the distal tip:

1st single blue marker: 75mm  
2nd single blue marker: 85mm  
3rd single blue marker: 95mm  
4th single blue marker: 105mm  
5th double blue marker: 120mm  
Mid-point of RED marker: 175mm  
Mid-point of wide BLUE marker: 285mm

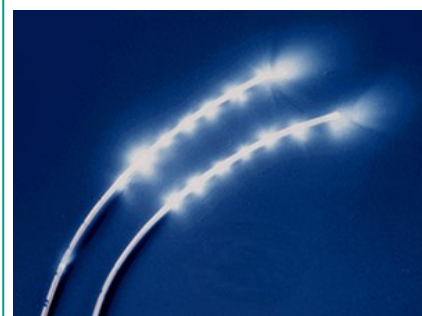

Description	Product code
<p><b>URIGLOW® Transilluminating Ureteric Stents</b></p> <p>Packed: Two stents in protective tray, complete with cystoscope bung and instruction sheet.</p> <p>For single use.</p>	R57412
<p><b>URIGLOW® Light Guide Coupler</b></p> <p>For attachment of URIGLOW® stents to fibre light sources with outputs &gt;250W.</p> <p>Storz fitting fibre cable. Reusable.</p>	R57411

For additional information, or to order one of our products, please call us on 01844 340 620 or email [info@aprmedtech.com](mailto:info@aprmedtech.com)

[aprmedtech.com](http://aprmedtech.com)  
[@aprmedtech](https://twitter.com/aprmedtech)

APR Medtech Limited  
The Sanderum Centre, 30a Upper High Street, Thame, Oxfordshire, OX9 3EX  
t. +44 1844 340 620 f. +44 1844 340 621 e. [info@aprmedtech.com](mailto:info@aprmedtech.com)



Rocket Medical PLC  
Rocket Medical plc., Sedling Road,  
Washington Tyne & Wear, NE38 9BZ



**Rocket**medical