

GATE OESOPHAGEAL


SELF-EXPANDING NITINOL STENT

Indicated for the management of malignant and benign oesophageal strictures

Dumbbell design to minimise risk of stent migration

Hand-woven nitinol hook and cross stent design provides optimal radial and low axial force

Gold radiopaque markers provide excellent visualisation under fluoroscopy

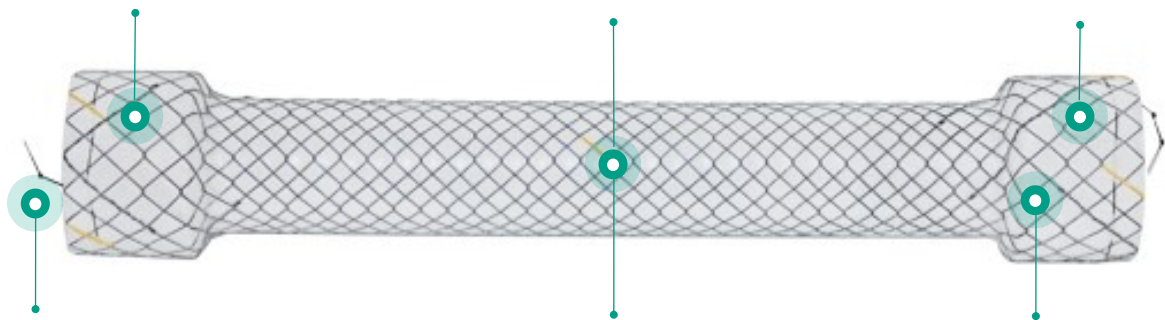

Retrieval string at both ends for both standard and "inside-out" removal


Results in minimal stent foreshortening upon deployment

Silicone coating designed to resist tissue ingrowth and facilitate safe removal


18Fr low profile delivery system aids placement without the need for predilation




Lumen conformity

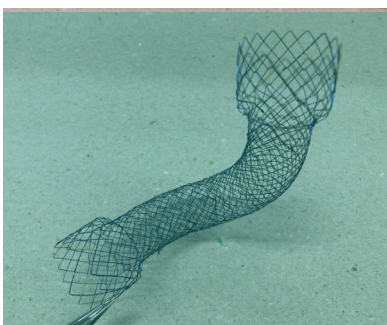

Radial strength

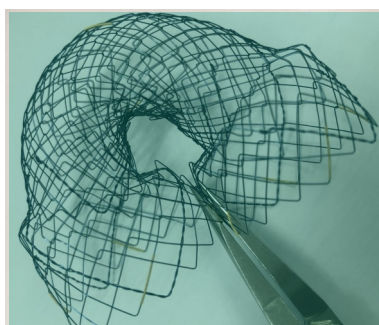

Full size range

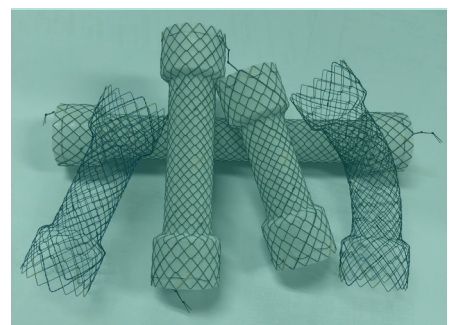


DESIGN CHOICE

The Gate Oesophageal Stent is available in a broad range of sizes. Please contact us with your preferred stent specification and we will provide you with additional ordering information.

STENT TYPE	DELIVERY SYSTEM		STENT	
	Diameter (Fr)	Length (mm)	Diameter (mm)	Total stent length (cm)*
UNCOVERED	18	700	18, 20, 22	8, 9, 10, 11, 12, 13, 14, 15, 16
COVERED	18	700	18, 20, 22	8, 9, 10, 11, 12, 13, 14, 15, 16

*Total stent length includes 2cm at each end to incorporate the flared / dumbbell design.

For additional information, please call us on [01844 340 620](tel:01844340620) or email info@aprmedtech.com

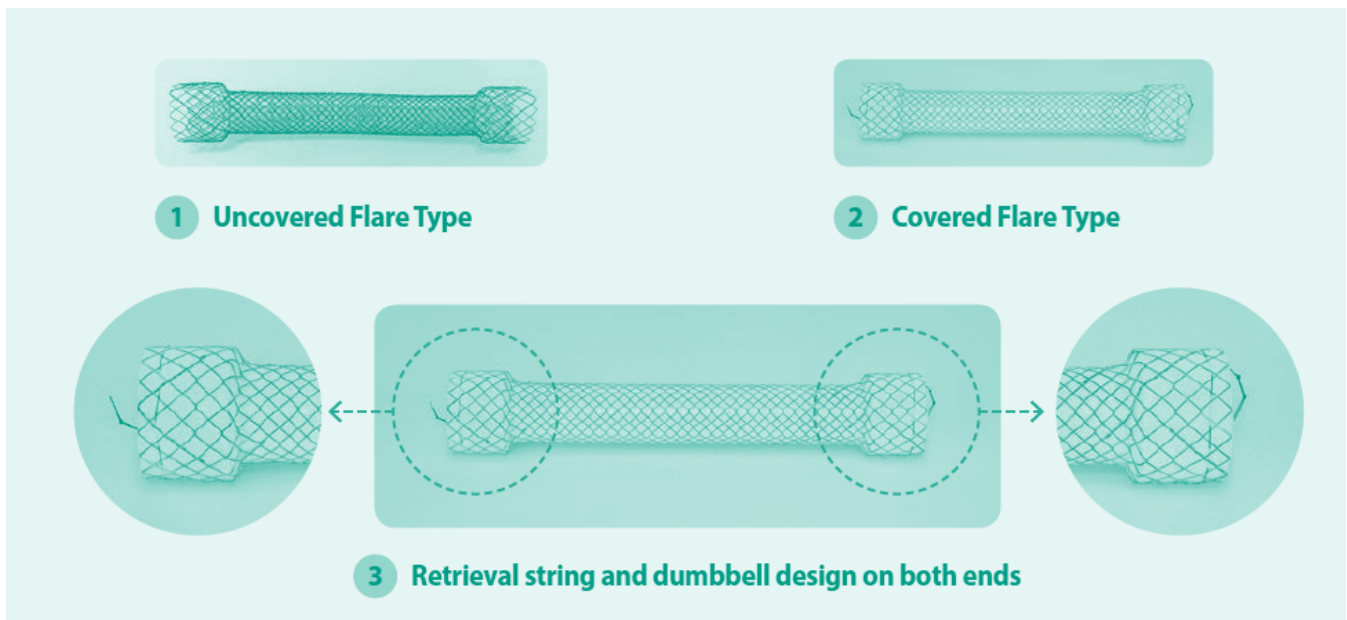

Manufactured by:



www.mediwood.co.kr

Supplied exclusively in the UK by:



aprmedtech.com

@aprmedtech

APR Medtech Limited

The Sanderum Centre, 30a Upper High Street, Thame, Oxfordshire, OX9 3EX

t. +44 1844 340 620 f. +44 1844 340 621 e. info@aprmedtech.com

