



GATE PYLORIC / DUODENAL

SELF-EXPANDING NITINOL STENT

Woven stent construction with small unfixed cells offers high flexibility with minimal shortening

The **nitinol weave** rotates in both planes and in different directions for optimal conformability


Gold radiopaque markers provide excellent visualisation under fluoroscopy

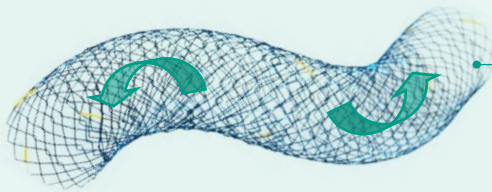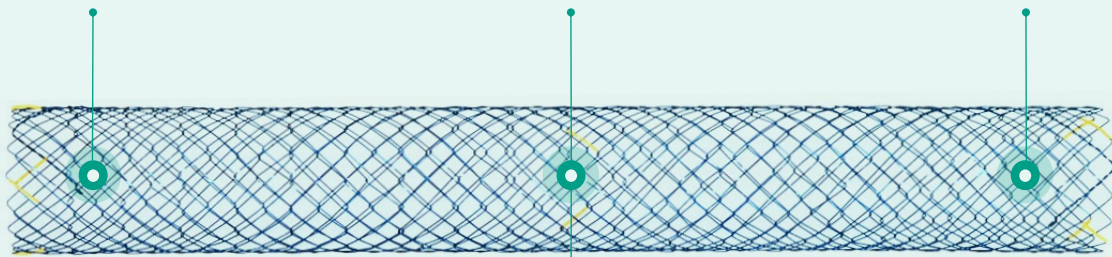

- Option of **flared or straight ends**
- Simple to use **10F TTS delivery system**
- Available in **starting diameter of 16mm**

Uncovered stent


Indicated for the management of malignant pyloric and duodenal strictures

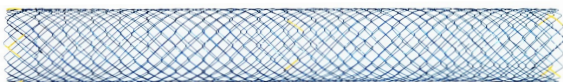

- **High flexibility** with minimal shortening
- **Excellent conformability and flexibility** to accommodate tortuous anatomy

Covered stent

Indicated for the management of malignant and benign pyloric and duodenal strictures

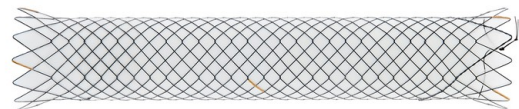

- **Silicone coating** designed to resist tissue ingrowth and facilitate safe removal
- **Retrieval string on distal end** to aid safe stent removal

DESIGN CHOICE

The Gate Pyloric / Duodenal Stent is available in a broad range of sizes, with three delivery system lengths. Please contact us with your preferred stent specification and we will provide you with additional ordering information.


STENT TYPE	DELIVERY SYSTEM		STENT	
	Diameter (Fr)	Length (mm)	Diameter (mm)	Length (cm)
UNCOVERED / STRAIGHT	10	1200, 1800, 2200	16, 18, 20, 22	4, 5, 6, 7, 8, 9, 10, 11, 12
COVERED / FLARED	10	1200, 1800, 2200	16, 18, 20	4, 5, 6, 7, 8, 9, 10, 11, 12

For additional information, please call us on [01844 340 620](tel:01844340620) or email info@aprmedtech.com


Lumen conformity



Shape flexibility

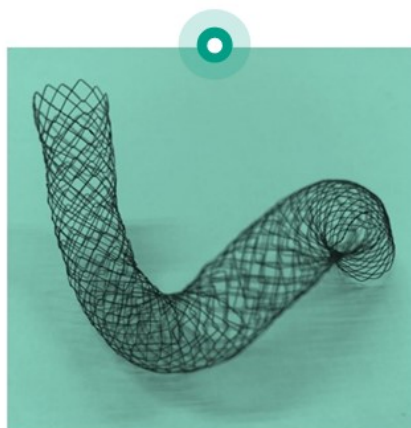

Radial strength

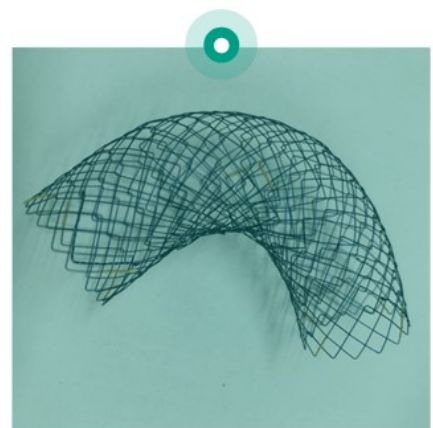

Manufactured by:



www.mediwood.co.kr

Supplied exclusively in the UK by:



aprmedtech.com
@aprmedtech

APR Medtech Limited
The Sanderum Centre, 30a Upper High Street, Thame, Oxfordshire, OX9 3EX
t. +44 1844 340 620 f. +44 1844 340 621 e. info@aprmedtech.com

