

## Interventional Cardiology Accessories



## Company overview

### OUR COMPANY

APR Medtech is a specialist independent medical technology company providing high quality medical devices and support services to the NHS and private healthcare sector.

Our products offer intelligent clinical solutions designed to help advance healthcare and improve patient outcomes. Our extensive experience within the healthcare industry, deep knowledge of the clinical areas in which we work and a genuine passion for what we do, ensures the delivery of real value to our customers.

### SERVICE AND SUPPORT

One of the key elements of our customer support service is product training. We want to be sure that our products are always used in the optimal way to the benefit of both our customers and their patients.

If you would like one of our team to visit, either to attend a procedure or run a training session on one of our products, then please get in touch. We will be more than happy to oblige.

### GET IN TOUCH

For additional information please call us on 01844 340 620 or email [info@aprmedtech.com](mailto:info@aprmedtech.com).

Additional information is also available on our website: [www.aprmedtech.com](http://www.aprmedtech.com).

Follow us on Twitter [@aprmedtech](https://twitter.com/aprmedtech) 



# VASCULAR ACCESS & INFLATION

## Arthesys Y-Connector


Arthesys Y-Connector (previously SuperKetch) is the original vascular access system designed for single-handed use.

**UNIQUE 'PUSH-PULL' DESIGN** for fast, accurate, single-handed device insertion

**SIGNIFICANTLY REDUCED BLOOD LOSS** provides a cleaner working environment and stable blood pressure monitoring

**ESSENTIAL FOR DES** offering safe insertion of drug-eluting stents and respecting the integrity of the coated surface

**EFFECTIVE FOR COMPLEX PROCEDURES** including bifurcation stenting, kissing balloons and multiple wire insertion

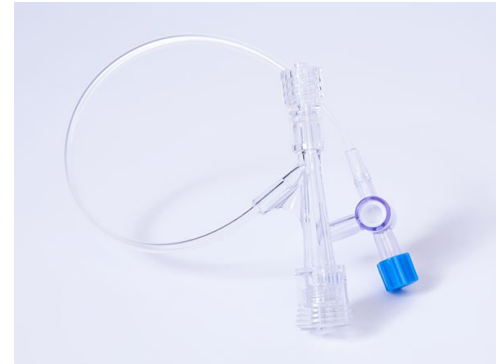

Manufactured by:  
**ARTHESYS**

### TECHNICAL SPECIFICATION

Component	Specification	Detail
Arthesys Y-Connector	Inner lumen compatibility	7F / 0.092" / 2.33mm maximum
Metallic insertion tool	External diameter	2.7F / 0.90mm
	Inner lumen diameter	1.9F / 0.025" / 0.64mm
	Usable length	10cm
Plastic peelable insertion tool	Length	9cm
	Inner diameter	3.4F / 0.044" / 1.14mm
Torquer	Guidewire compatibility	From 0.010" minimum to 0.035" maximum

### ORDERING INFORMATION

Product code	Description	Box qty
18020011 (KESU0002)	Arthesys Y-Connector with 20cm side arm. Guidewire introducer & torquer 0.010" / 0.014"	10
18020012 (KESU0005)	Arthesys Y-Connector with 50cm side arm. Guidewire introducer & torquer 0.010" / 0.014"	10
18020007 (KESU0004)	Arthesys Y-Connector with no side arm. Guidewire introducer & plastic peelable insertion tool.	10

## A125 Inflation Device

The A125 inflation device combines robust ergonomic design with high build quality to deliver excellent operator control and reliable, consistent performance.

**CONTROLLED INFLATION** with quick release mechanism

**LUMINESCENT DIAL** for 'lights down' interventional procedures

**30ATM** pressure gauge

### ORDERING INFORMATION

Product code	Description	Box qty
A125	30atm Inflation Device	5



Manufactured by:  
**KDL**

# THROMBUS ASPIRATION

## *Pegase HYDRO* Thrombus Aspiration Catheter

Keep it simple and efficient

**CONSISTENT SHAFT CONSTRUCTION** and removable stylet to avoid choosing between crossability and kink-resistance

**TRANSLUCENT** aspiration lumen for an easy control after catheter withdrawal

**SINGLE-HANDED** operator device

**SINGLE RADIOPAQUE MARKER** mounted on a short 0.015" distal tip for accurate positioning

**HIGH RATIO ID/OD** enabling outstanding aspiration speed



### TECHNICAL SPECIFICATION

Coating	Hydrocoat® on distal part
Usable catheter length	145cm
Entry profile	0.015"
Extraction lumen	0.93mm <sup>2</sup> for 6F / 1.43mm <sup>2</sup> for 7F
Aspiration speed	12.5s for 6F / 7s for 7F
Radiopaque marker	1mm at 4mm from the tip
Guidewire compatibility	0.014"



### ORDERING INFORMATION

Product code	Description	Box qty
04030105	Pegase HYDRO Aspiration Catheter 6F compatible	1
04030106	Pegase HYDRO Aspiration Catheter 7F compatible	1

Manufactured by:



For additional information, or to order one of our products, please call us on 01844 340 620 or email [info@aprmedtech.com](mailto:info@aprmedtech.com).

Proud to work in partnership with:



[www.arthesys.com](http://www.arthesys.com)



[www.new.kdlchina.com](http://www.new.kdlchina.com)

Distributed in the United Kingdom by:



[aprmedtech.com](http://aprmedtech.com)  
[@aprmedtech](https://twitter.com/aprmedtech)

APR Medtech Limited  
The Sanderum Centre, 30a Upper High Street, Thame, Oxfordshire, OX9 3EX  
t. +44 1844 340 620 f. +44 1844 340 621 e. [info@aprmedtech.com](mailto:info@aprmedtech.com)

